SMART BATHROOM SCALE

Model:

HBS-BTH180WH

- Bathroom scale
- Capacity: 180 kg
- Bluetooth function

HEINNER

1. INTRODUCTION

Please read the instructions carefully and keep the manual for future information.

This manual is designed to give all necessary instructions concerning installation, use and maintenance of the unit. In order to operate the unit correctly and safely, please read this instruction manual carefully before installation and use.

2. CONTENT OF YOUR PACKAGE


- ➔ Bathroom scale
- ➔ User manual
- ➔ Certificate of warranty
- ➔ Declaration of conformity

3. FEATURES

GENERAL


- 1.To ensure the accuracy of the data, it is recommended that each measurement be carried out at the same time of day. Don't use in strong electro-magnetic condition and field with intense shake.
- 2.Stand on the scale steadily and touch the electrodes on the platform with your bare feet. Tip over will lead to a fall when you step on the edge of one side of the scale. Keep standing firmly on scale when in use.
- 3.Do not use the scale in humid, overheated environment. (Temperature range from 5°C to 35°C)
- 4.When the battery runs out, please remove the battery in time to prevent leakage of the battery lead scale body corrosion.
- 5.If the product cannot work properly, then - For the scale using 1.5V AAA batteries, please check whether battery is installed, or battery power is low. Inset or replace new batteries. - For the scale with built-in rechargeable battery, please press the restart button to reset the scale with the aid of a sharp object.
- 6.Do not drop, shock, or strike the scale as it is a precise instrument.
- 7.Keep the scale clean. Do not use acid cleaning agent.

DISPLAY INDICATIONS

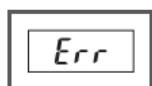
OVERLOAD INDICATION

The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage.

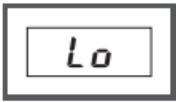LOW BATTERY INDICATION

Battery power is low. Please install new battery

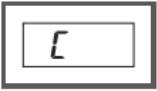
CALIBRATION

Display "C" when you step off the scale, means this weight reading with error. Please re-test

4. USING THE APPLIANCE

Please set the weight unit by APP. If unavailable, you can press UNIT button on back of scale to convert the weight unit.

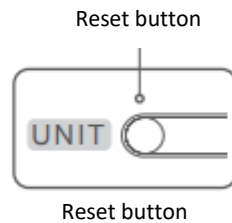Remark:

- 1.The unit system of scale depends on scale destination or the requirement of distributor.
- 2.Product with rechargeable battery is with a restart button as shown in the picture. If the product cannot work properly, please press the restart button to reset the scale with the aid of a sharp object.

DOWNLOAD AND SET UP THE APPLICATION

Search "MovingLife" in "App Store" or "Google Play" to download the Application. Or in another most convenient way by scanning the QR code below.


App: MovingLife

Tips : App is keeping updating, so it may differ from showed on manual. For the final and accurate information, please check on App.


REGISTERING AND CONNECTING

Open App "MovingLife", connect the device with App and submit your personal profile. (Before connection, please make sure Bluetooth is turned on and scale LCD is lighted up) Note: Please make sure your set personal information is correct, in order to get accurate measurements.

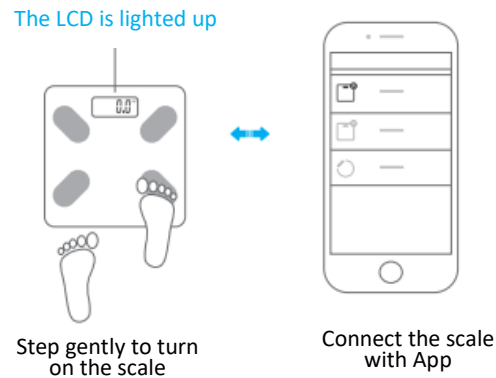START WEIGHING

Place the scale on the hard and flat floor. Stand on the scale with bare feet. The weight reading is stable and flashes, then begins to be transmitted. The scale shuts off after several seconds. App will show and save all indexes and history record.

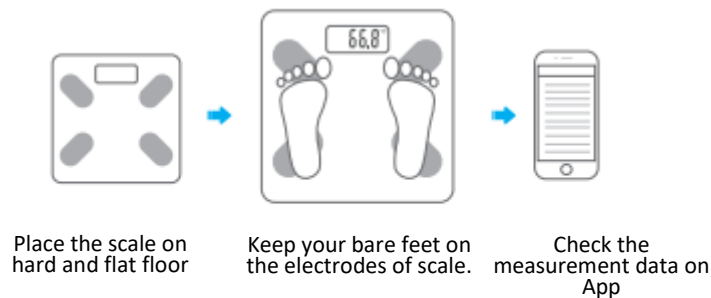Note: To make a proper measurement of body fat, please keep your bare feet on the electrodes of scale.

ABOUT THE BODY FAT SCALE


The product is a smart body fat scale with a high-precision strain gauge sensor, which can help you to store and keep track of your weight on smart phone. This scale can transmit your weight data to your smart phone via Bluetooth technology, as well as make graph of your weight results, from which you are able to keep track of the trend anytime. Download the App on your phone to use Smart scale, you will start your first step of weight control.

USER-FRIENDLY TIPS

The product can be used to weigh and calculate body fat % & hydration % after you set gender/ height/age parameters etc. These measurements are only indicative of body fat status. This product is not for medical treatment and can't be used in any medical institution. Growing children, pregnant women, those with symptom of edema and those on heart bypass treatment are not suggested to use the product.

MULTIPLE USERS SUPPORTED

The scale can be connected with multiple users.

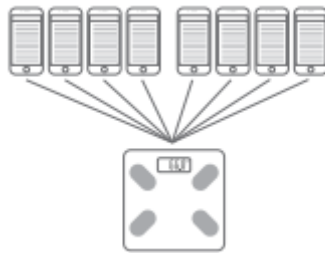Note: Final user amount depends on the product


THE DISPLAY OF AMBIENT TEMPERATURE

This product can detect and display ambient temperature automatically if starting up. And it would constantly refresh the latest ambient temperature until the power off. If you change the service environment of the product, it is recommended to check the temperature display 30 minutes later to ensure its accuracy. (Some products don't have temperature display, please refer to the actual one.)

5. BATTERY

If it's powered by AAA batteries, please install as shown below: First Use: Open the battery cover on the back of scale and put out the insulating strip.

Change battery: Open the battery cover on the back of scale, then install the batteries properly according to the pictures shown below.

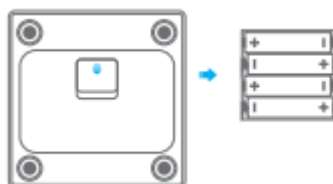Note : After battery installation, please close the battery cover. Please use the corresponding battery. Remove battery if it is not used for a long time (more than three months).

If the scale is with built-in rechargeable battery

1. The scale should be fully charged when firstly used or after long-time stored. Battery symbol flashes during charging and displays constantly when charging is finished. It takes about 2 hours to finish.
2. Please use android 5V charger with current equal to or smaller than 2A. DO NOT use a charger with voltage larger than 5V.
3. Suggest charging every 3 months if the scale is not used for a long time.

Remark: Battery symbol is shown on LCD during charging. For your safety, please do not use the scale during charging. Only start weighing after charging is finished and removing charger.

6. TROUBLESHOOTING

The software can't obtain weighing results from Smart scale.

- Please confirm the connection sign on the App interface to know the connection status.
- Confirm the Bluetooth is on in your device and the Smart scale is on within the valid range of signal coverage.
- Check whether the battery power is running low.
- Relaunch the application.
- Turn off the Bluetooth and on again, then bridge the Bluetooth connection between scale and smart phone again.

Fail to connect to your device.

- Please make sure only one device is setting up connection with the Smart scale.
- Please make sure your device is equipped with Bluetooth 4.0 or above.

LCD has no response.

- If it's powered by AAA batteries, please first check the batteries, and replace them if necessary. If the scale is with built-in rechargeable battery, please use the USB cable to recharge the scale if necessary.
- Repairs may only be carried out by our authorized customer service or dealers.

CUSTOMER CARE AND SERVICE

Always use original spare parts.

When contacting our Authorised Service Centre, ensure that you have the following data available: Model and Serial Number.

The information can be found on the rating plate. Subject to change without notice.

TECHNICAL DATA

Capacity	180 kg
Power	Battery type: 3x1.5V AAA (not included)

Thank you for purchasing this product. If you need support with your product, visit our website using the links below.

Get User manuals: <https://www.heinner.ro>

Get Service information: <https://www.heinner.ro>

**Environment friendly disposal**

You can help protect the environment!

Please remember to respect the local regulations: hand in the non-working electrical equipment's to an appropriate waste disposal center.

HEINNER is a registered trademark of **Network One Distribution SRL**. Other brands and product names are trademarks or registered trademarks of their respective holders.

No part of the specifications may be reproduced in any form or by any means or used to make any derivative such as translation, transformation, or adaptation without permission from NETWORK ONE DISTRIBUTION.

Copyright © 2013 Network One Distribution. All rights reserved.

www.heinner.com, <http://www.nod.ro>

This product is in conformity with norms and standards of European Community.

Importer: **Network One Distribution**

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Copyright © 2013 Network One Distribution. Toate drepturile rezervate.
Copyright © 2013 Network One Distribution. Всички права запазени.
Copyright © 2013 Network One Distribution. Minden jog fenntartva.